

# FD-PHARMA™ SERIES II DIE CUT


## Liquid Unit Dose Cup Packaging System

Manufactured by Medical Packaging Inc.

Liquid Unit Dose Cup  
Packaging Made Simple!



Medical Packaging Inc. (MPI)'s **FD-Pharma™ Series II Die Cut**, a dual-head liquid unit dose cup packaging system that prints, seals and cuts two (2) cups at a time, represents a low cost packaging solution to small pharmaceutical manufacturers and CDMO's seeking to expand their market base with new and expanded packaging capabilities. **FD-Pharma™ Series II Die Cut** can quickly, easily, and inexpensively add liquid cup packaging capability to your manufacturing operations. For the CDMO, adding the **FD-Pharma™ Series II Die Cut** can provide additional value to your current clients who are seeking to minimize the number of packaging partners they are engaged with as well as attract new clients by offering liquid cup packaging capability not previously available.

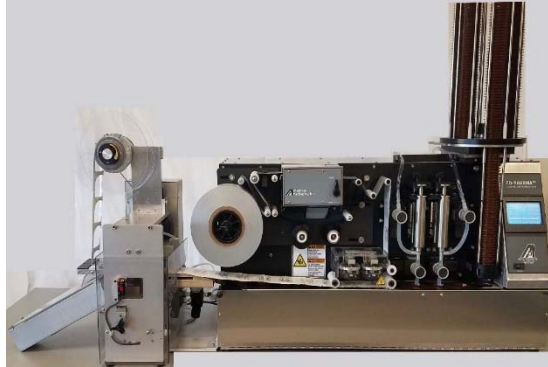

MPI's **FD-Pharma™ Series II Die Cut** was developed to address the void that exists for early liquid product development and small batch liquid validation runs as well as lower volume production needs not suited for high speed packaging lines. **FD-Pharma™ Series II Die Cut** provides the ability to support early phase development, small batch runs, and medium production runs. **FD-Pharma™ Series II Die Cut** offers the opportunity to validate a new liquid drug using production equipment that is offered at a much lower cost production machine. Once the liquid drug product has gained significant market share that warrants higher speed production capabilities, production can easily be transitioned to the higher speed packaging line.

MPI's **FD-Pharma™ Series II Die Cut** offers the perfect packaging solution for lower volume products that represent good margins but are too low in volume to be run profitably on the higher speed packaging lines used for blockbuster products. **FD-Pharma™ Series II Die Cut**, with its quick changeover and easy to operate features, would allow for multiple small to medium profitable production runs of varying liquid products on one machine.

*MPI's **FD-Pharma™ Series II Die Cut** is the only liquid cup unit dose packaging machine available that is able to produce FDA-compliant packages, addressing initial product development to medium production volumes, at a low price point, producing high quality, professionally packaged unit dose cups rivaling those produced on much higher cost packaging equipment!*

*\*The **FD-Pharma™ Series II** system can be retrofitted with the die cut module*

### KEY FEATURES

- *Small footprint, small price point, with professional FDA-compliant results for precise liquid dosing in a convenient unit dose cup.*
- *The **FD-Pharma™ Series II Die Cut** is well suited for clinical trials, small production, to full production mode . . . all in one machine!*
- *Easy to operate by line technicians with minimal training required to be proficient.*
- *Quick changeover with minimal cleaning required.*
- *Cost of ownership is very low with the **FD-Pharma™** due to very low maintenance costs for operating and maintaining the machine.*
- *Produces professional results at a low cost per dose (CPD).*
- *Barcoding and serialization capabilities (when using **Pak-EDGE™ UD**) are standard features.*
- ***FD-Pharma™ Series II Die Cut** with optional scanner provides 100% verification of a readable barcode on the unit dose package.*
- *Pumping System: Dual Stainless Steel Syringe, positive displacement, ideal for high viscosity fluids.*
- *Standard cup sizes available: 7mL, 15mL, 25 mL 35mL*

### TECHNICAL SPECIFICATIONS

Dimensions:	Package: 37"L x 14"W x 18"H Die Cut: 22"L X 12"W x 24"H
Weight:	Package: 130 lbs. Die Cut: 40 lbs.
Power:	Package: 115VAC, 50/60Hz, 500VA Die Cut: 115VAC 50/60Hz, 50VA
Sealing Temp:	220° ± 5° C
Rate:	Up to 23 packages per minute (ppm)
Fill Accuracy:	+/- 0.5%



To request information, contact:  
[pharmasales@medpak.com](mailto:pharmasales@medpak.com)

FD-Pharma and Pak-EDGE  
are trademarks of Medical Packaging Inc.

**Medical  
Packaging Inc.**

8 Kings Court | Flemington | New Jersey | 08822 | USA  
Toll-Free.800-257-5282 | Phone.609.466.8991 | Fax.609.466.3775 | [www.medpak.com](http://www.medpak.com)



050520

# FD-PHARMA™ SERIES II DIE CUT

## Liquid Unit Dose Packaging System

Manufactured by Medical Packaging Inc.

# OPTIONS



### Barcode Verifier Option

SCANS FINAL UNIT DOSE PACKAGE  
TO VERIFY READABLE BARCODE



### Air Rinser

TARGETS A SHORT BURST OF AIR INTO CUP  
DIRECTING DEBRIS INTO  
COLLECTION NETTING (NOT SHOWN)



### Auto-cup Indexer

HOLDS AND FEEDS UP TO 300 SMALL CUPS  
OR 180 LARGE CUPS

### Color Cups

CHOOSE FROM STANDARD OR SPECIAL  
ORDER CUSTOM COLORS



CE

To request information, contact:  
[pharmasales@medpak.com](mailto:pharmasales@medpak.com)

 **Medical  
Packaging Inc.**

8 Kings Court | Flemington | New Jersey | 08822 | USA  
Toll-Free.800-257-5282 | Phone.609.466.8991 | Fax.609.466.3775 | [www.medpak.com](http://www.medpak.com)



050520

FD-Pharma and Pak-EDGE  
are trademarks of Medical Packaging Inc.