

AUTO-WRAP[®]

Syringe Labeling System

Manufactured by Medical Packaging Inc., LLC



DECREASE PACKAGING COSTS

INCREASE OPERATIONAL EFFICIENCY

REDUCE MEDICATION ERRORS

The Medical Packaging Inc., LLC (MPI) **Auto-Wrap[®] Syringe Labeling System** is an easy-to-use solution that automates the process of labeling syringes. This semi-automatic system produces a printed label and automatically wraps the printed label around the syringe in one single motion.



MPI's Auto-Wrap[®] Syringe Labeling System can accommodate any type of oral slip or injectable syringe from ½ mL up to 60 mL without any need for adjustment or changing of parts. The open face design and easy-to-use controls allow for easy loading of label stock and changeover.

Flagged and tacked labels are offered with print areas ranging from 1" – 1.5". Color flag labels can be used for various purposes, such as identifying pediatrics, variants/doses of a drug, storage conditions etc.

KEY FEATURES



- TIME SAVINGS** - An independent time motion study conducted on MPI's Auto-Wrap[®] concluded that time savings of 1 ½ - 2 ⅓ minutes per batch of 25 syringes can be achieved by implementing MPI's Auto-Wrap[®] Labeling System.

Key Metrics	Syringe Size	Time Savings (per 25 Syringes)
Straight labeling	Any	1 1/2 minutes
Folded labeling	1 mL	2 2/3 minutes
Folded labeling	> 5 mL	2 1/3 minutes

Study conducted by Summerfield Consulting


- COST SAVINGS** - Low cost per dose (CPD) . . . ~ 4 cents/label.
- PRODUCTIVITY** - Increase production throughput using automation while eliminating manual wrapping which can contribute to repetitive stress injuries (RSI). *To assure maximum adhesion and legibility, MPI's Auto-Wrap[®] Syringe Labeling System requires the use of label media developed by MPI in conjunction with MPI's Pak-EDGE[®] UD Barcode Labeling Software 2.0, which is included with the labeling system.*

FLAGGED LABELS

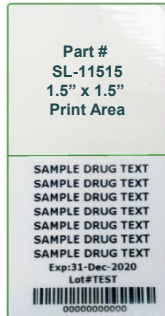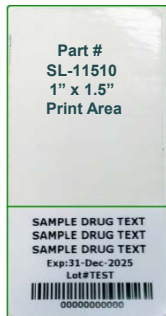

The top clear portion of the label is wrapped around the syringe with semi-permanent adhesive. The white/printed portion extends and hangs from the syringe and is deadened to reduce the tackiness.

TACKED LABEL

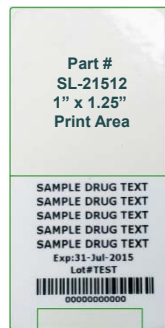

The top clear portion of the label is wrapped around the syringe with semi-permanent adhesive. The white/printed portion has a clear adhesive area which will keep the label tacked to the syringe. This will allow for better storage and for the label to be untacked from the syringe for viewing.

TECHNICAL SPECIFICATIONS

Dimensions:	17.5" W x 7.5" D x 11.0" H
Weight:	25 lbs.
Safety:	CE, NRTL (MET) US, C
Power:	115/230V, 50/60 Hz, 70 Watt

The **Auto-Wrap[®] Syringe Labeling System** is the perfect complement to **MPI's Auto-Draw[®] Syringe Filling System!**



To request information, contact:
sales@medpak.com



8 Kings Court | Flemington | New Jersey | 08822 | USA
Toll-Free.800-257-5282 | Phone.609.466.8991 | Fax.609.466.3775 | www.medpak.com



Auto-Wrap and Pak-EDGE are registered trademarks of Medical Packaging Inc., LLC.

101222